

PECS CAN INTERFACE WITH:-

- MATLAB
- LABVIEW
- MS Word
- MS Excel
- MS Access
- C++, net, C#,J++
- Visual Basic
- Java

Vicon Motion Systems has been producing motion capture and analysis systems since 1984. Vicon PECS allows you to streamline external program workflow within the Vicon Pipeline environment.

ABOUT VICON PECS


Vicon's latest Plug-in, PECS, turns Vicon Workstation into a Windows Active X Automation Server. This gives any application supporting ActiveX technology access to Workstation data and processing API.


Using the PECS Plug-in to Workstation, a wide variety of client applications, scripts and processes can be evaluated in the Vicon Pipeline. The Plug-In offers the additional capability to provide ActiveX access for processing of 3D marker data and analogue data within other applications.

PECS gives access to all trial data; including 3D reconstructions, analogue and force plate data, events and contexts; subjects and subject parameters.

PURCHASE VICON PECS ONLINE!

Please visit www.vicon.com and follow the links to purchase now by credit card. Alternatively please contact your local Vicon office or distributor.

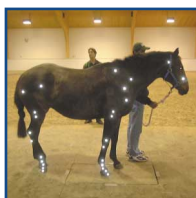

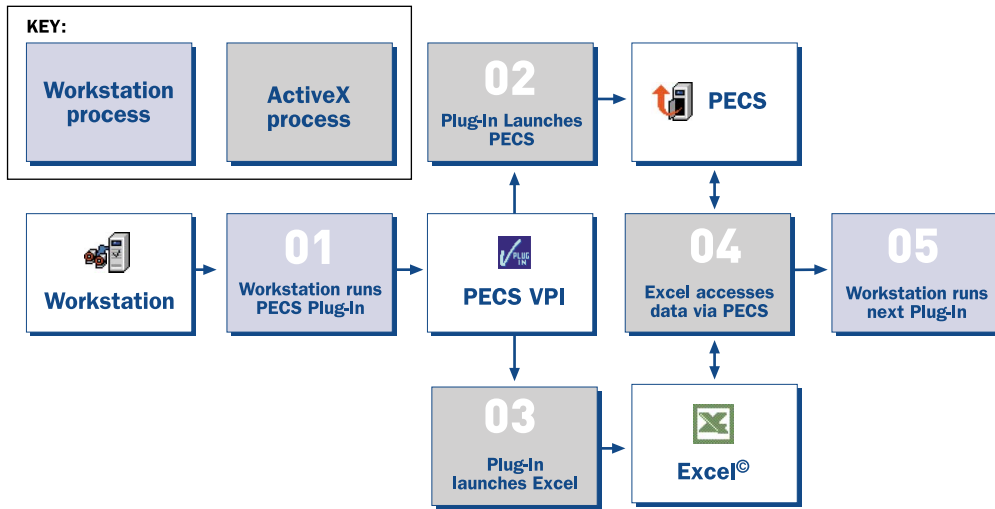

HIGHLIGHTS

Faster processing

PECS enables you to process externally written code without the need to export from the Vicon suite.


Plug-in Architecture

Since PECS is installed inside the Vicon Pipeline, your workflow can be accelerated by the use of Trial Types and Batch Processing.

SOME EXAMPLES

- Run your existing MATLAB® scripts in the processing pipeline
- Run Virtual Instruments on trial data using LabView™
- Automate reports using Word templates and Visual Basic
- Automate reports and calculate statistics in trial data using Excel
- Write Plug-in executables in other languages

ALL WITHIN THE VICON PIPELINE!

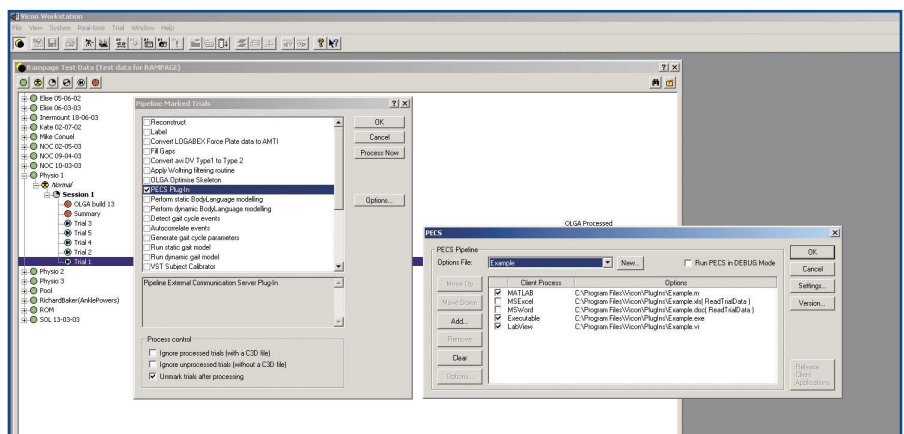

Vicon Motion Systems Inc - USA
 9 Spectrum Pointe Drive
 Lake Forest, CA 92630, USA
 T. +1 949 472 9140 F. +1 949 472 9136
 E. sales@vicon.com

Vicon Motion Systems Ltd - UK
 14 Minns Business Park, West Way
 Oxford, OX2 0JB, UK
 T. +44 (0) 1865 261800 F. +44 (0) 1865 240527
 E. sales@vicon.com